


TAP MILL

Dust contained milling system

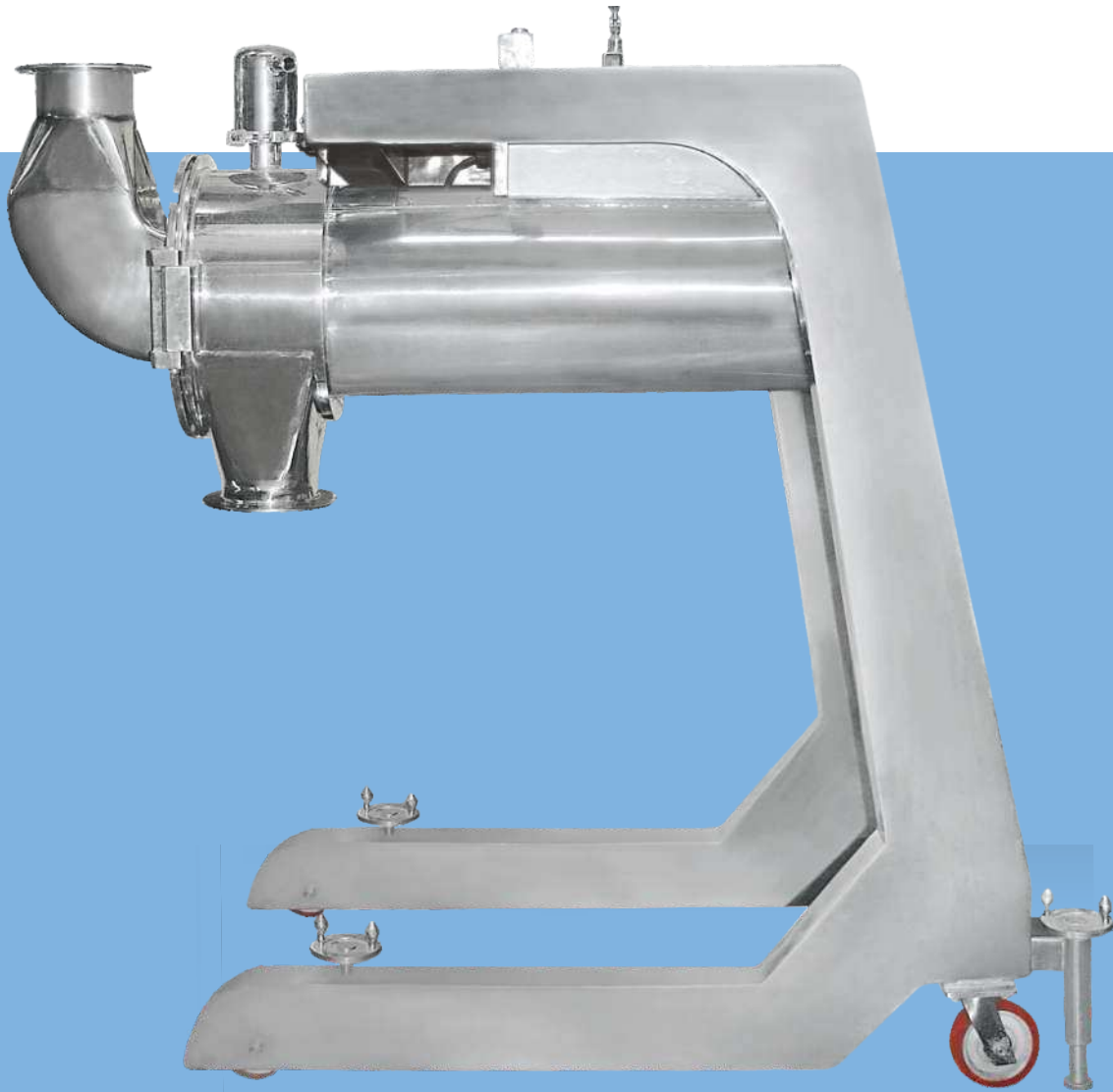

CE

TAPASYA

TAP MILL is used for high speed milling & sizing of dried granules, powder and a wide range of wet and dry materials without special attachments. This machine operates on the principal of variable force of swing cutting knives edges blades. Exclusively used in Pharmaceutical, Food, Chemical & Bio-technology Industries. As compared to the four common principles of size reduction which often do not produce controlled size reduction, this machine utilises the principle pf variable force swing hammer blades having both knife and impact edges rotating with a carefully selected screen to control size reduction.

TAP MILL

Dust contained milling system

APPLICATIONS

Milling pulverising, chopping, shredding.


Online dust contained milling.

TECHNICAL SPECIFICATION

Model	Working Cap in kgs/hr	RPM (NOMINAL)	HP
SAI-TM-1	20 - 50*	2880	1
SAI-TM-2	50-100*	2880	2
SAI-TM-3	200 - 500*	2880	5
SAI-TM-4	500 - 750*	2880	5

SALIENT FEATURES

- ! Easy to clean
- ! Minimum maintenance
- ! Bin charging system for Loading and Unloading of Powders and Granulesfor Complete Dust free operation
- ! The equipment is made from SS 316
- ! Excellent for Contained milling
- ! Capacity available from 50-2000 kgs/hr (Higher capacities upto 6000 kgs/hr also provided on demand)

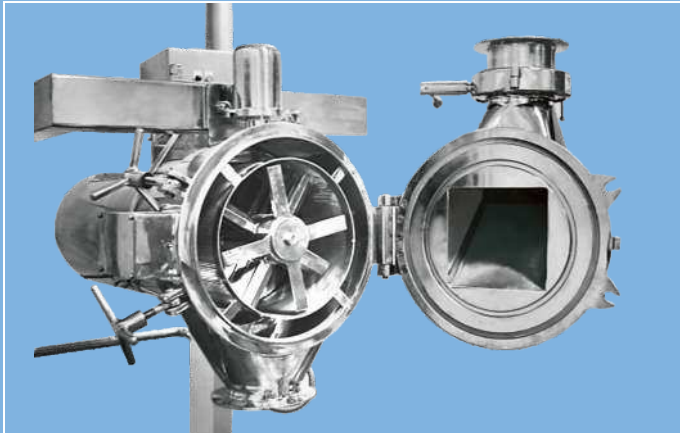

Specification and technical data are subject to change without prior notice



Manufactured & Exported By:

Tapasya Engineering Works Pvt. Ltd.

A/212, Road No.30, Wagle Industrial Estate, Thane (W) – 400 604. India

Tel.: +91 - 22- 6157 9400 / 2582 2287 • Fax : +91 - 22- 2582 5243

E-mail: sales@tapasyaindia.net • info@tapasyaindia.net • Web : www.tapasyaindia.net

Our Representatives in USA:

Daiichi Jitsugyo (America), Inc.

939 AEC Drive, Wood Dale, IL 60191

Phone : 630-776-2459 • Email: Rnelson@djausa.com

Web : www.dja-pharma.com • www.tapasyaglobalusa.com

